RUD-TJH Experiences

Jalil Hosseini

Professor of Urology Reconstructive Urology Dept. Tajrish Hospital, Tehran – Iran

President, A-P SRUAS
Past President, UAA

Honorary Member, EAU Past President, IUA











- History:
 - 90 years Old,
 - Long anterior U.Stricture
 - Several Dilatation, Comorbid
 - Every 3 month in Operation Room with Anesthesia


History:

- 85 Years Old
- Prostatectomy(or TURP)
- Membranous Stricture
- Several Dilatation
- Every 2 month in Operation Room with Anesthesia









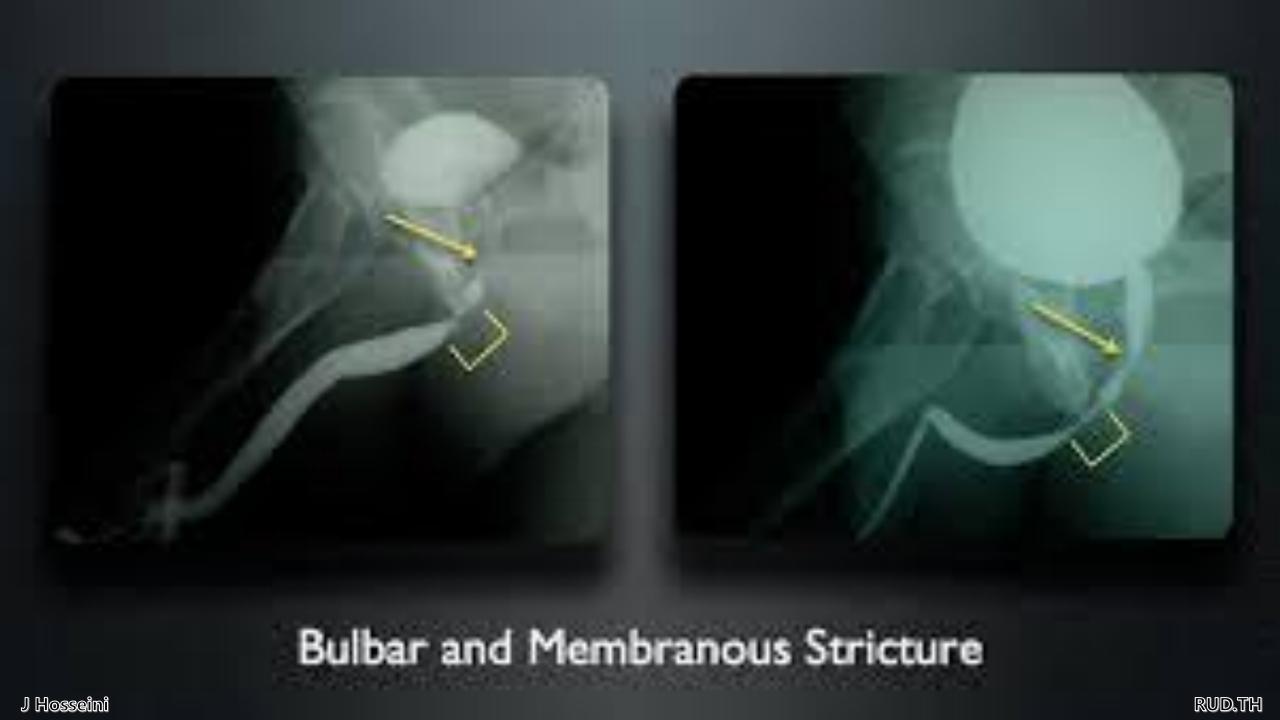






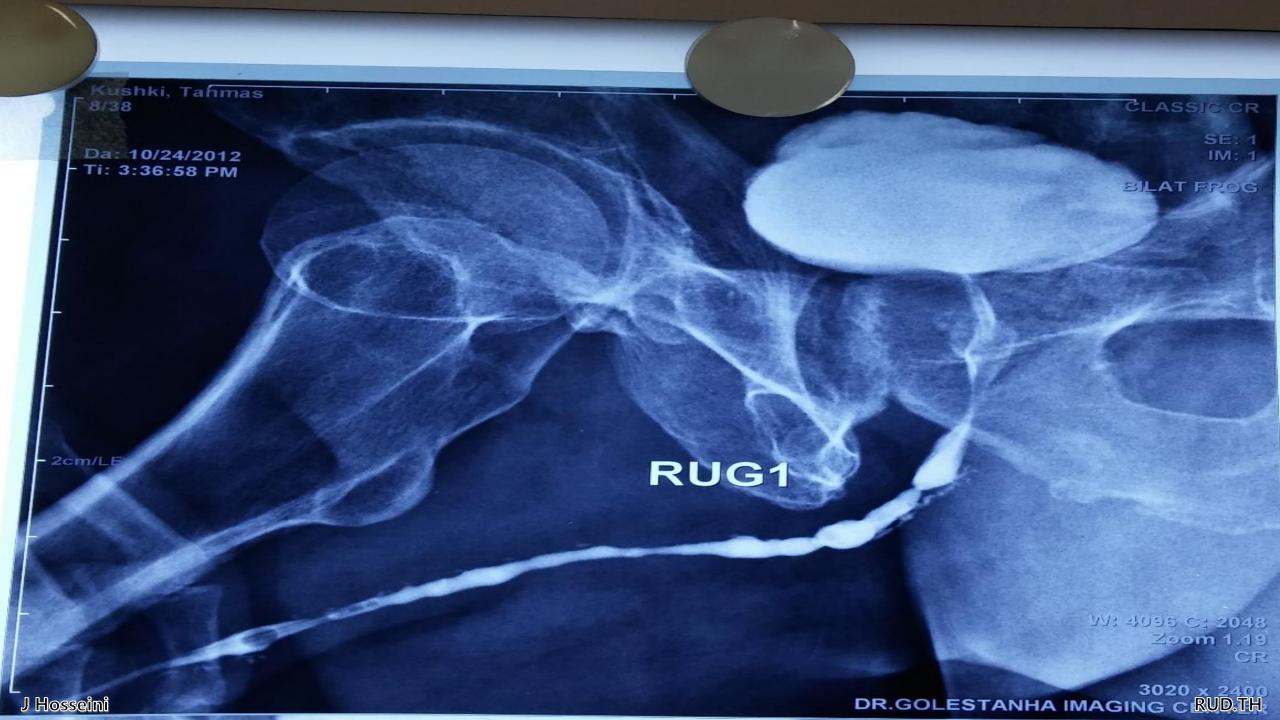

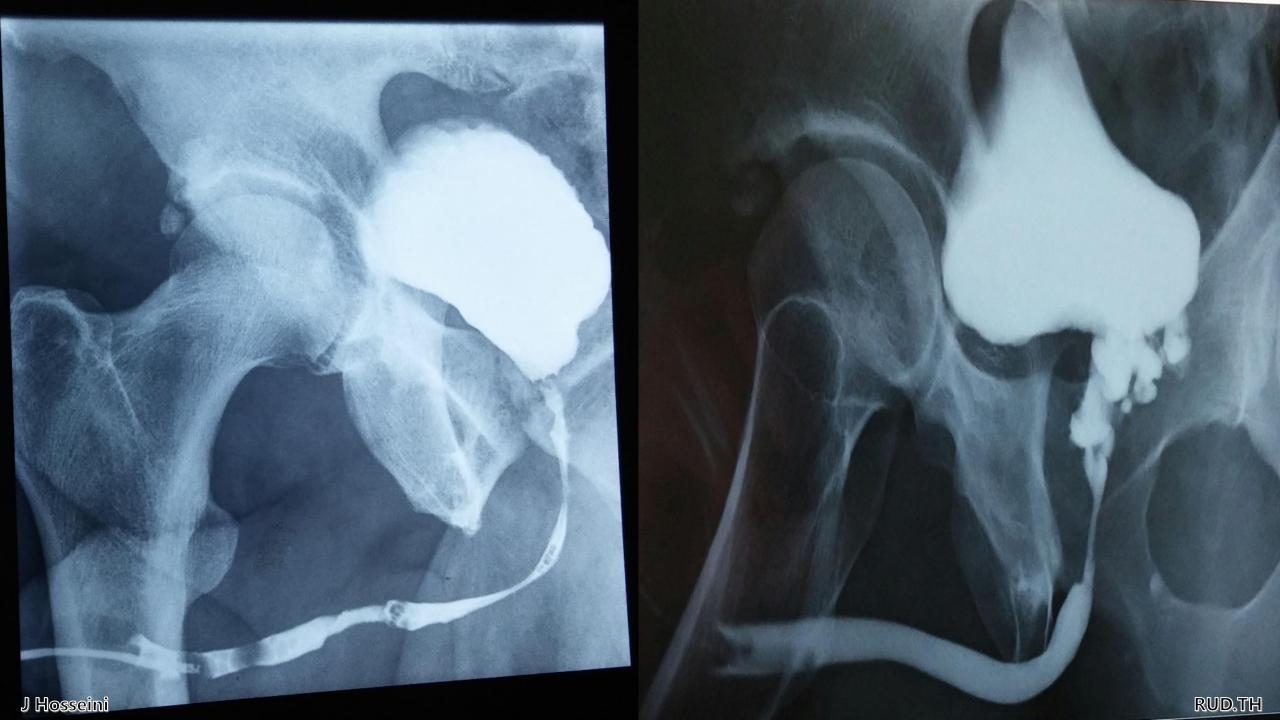

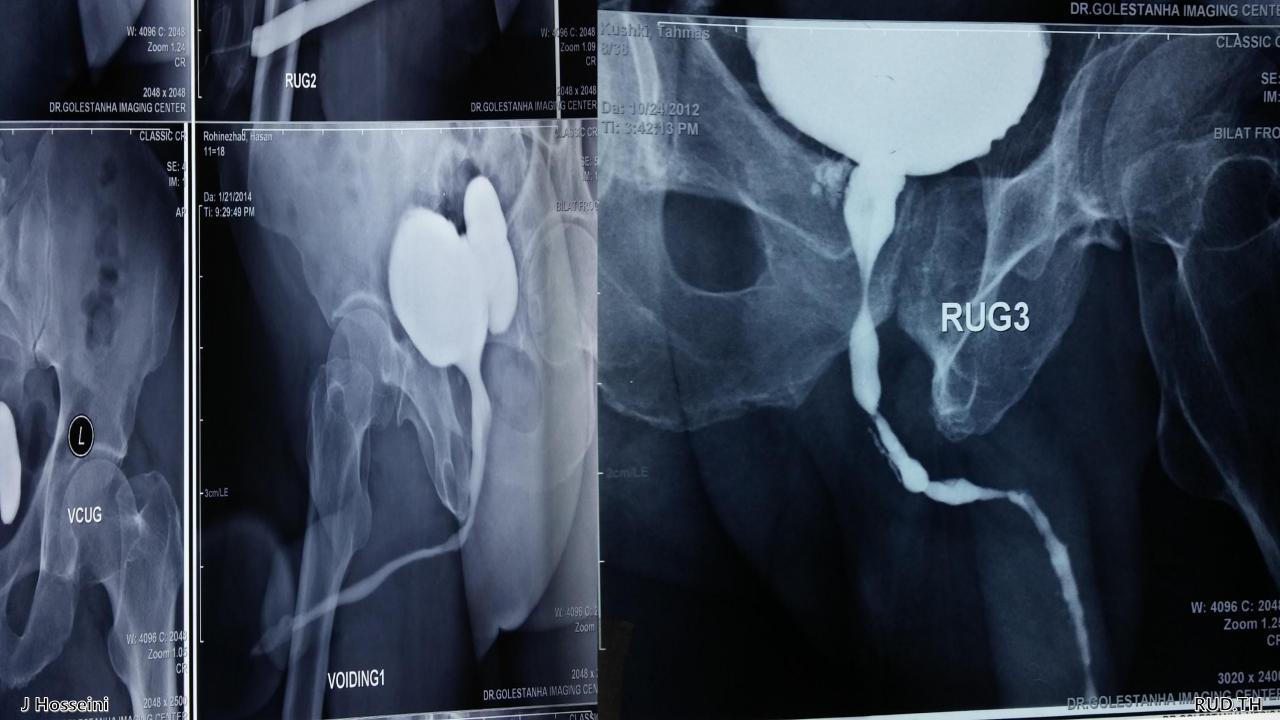


- History:
 - 85 Years old 90 Years Old
 - Many operation room admit

-How many time do you want to repeat

Hydraulic Urethral Dilatation

- One of the effective therapies for urethral stenosis is the use of a urethral hydraulic catheter.
- In order to take advantage of the optimal and appropriate therapeutic response
- With proper use of this catheter, the patient will not require frequent intraoperative delivery, anesthetic tolerance and recurrent dilation.
- The patients need Primary urethral dilatation to accommodate HUD
- Primary education is done at the clinic in the presence of the urologist

HUD catheter

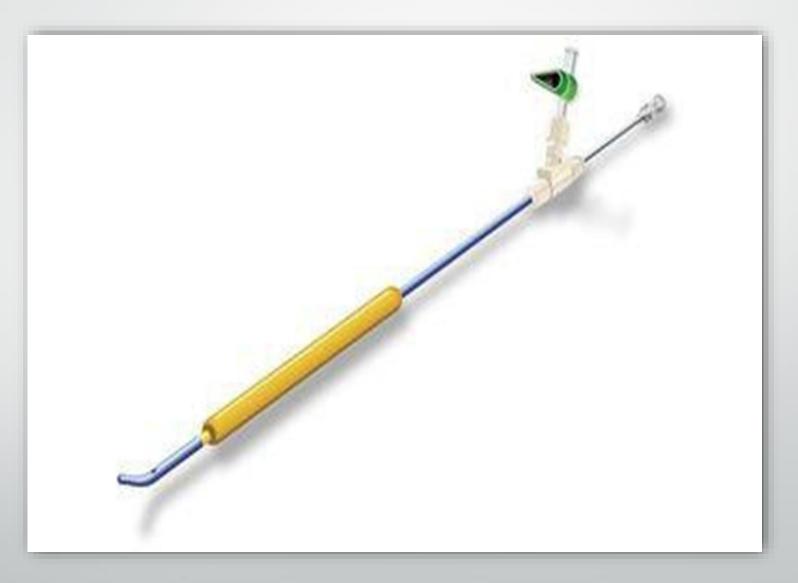

Technical Points

First, remove the catheter from the compartment and clean it with a sterile serum.

• Before inserting the catheter, it is important to use the lubricants

- The catheter should be inserted into the urethra with its tilting upwards and note that the entire yellow part of the catheter is inside the urethra.
- It is recommended to use distilled water to fill the balloon around the catheter and fill the balloon to the extent that the patient feels pain (4 to 8 cc).
- Depending on the urologist's advice, the duration of the catheter in the urethra will be between 4 and 6 hours.
- After using the catheter, the catheter is thoroughly washed with a clean wash, sterilized in a sterile cleanser and kept in a container in which a single formalin tablet is inserted.

•

• At the time of re-use with the prescribed urologist, the steps above should be repeated from the beginning.

